


Sample Pilot Purge Gas Filter Assembly

This filter connects to the Sample Pilot and provides the sterile barrier at the purge gas connection. It must be replaced every time a bioreactor has finished a campaign. The purge filter is assembled to the sample pilot prior to autoclaving. The purge filter helps ensure contamination-free operation of the bioreactor system.

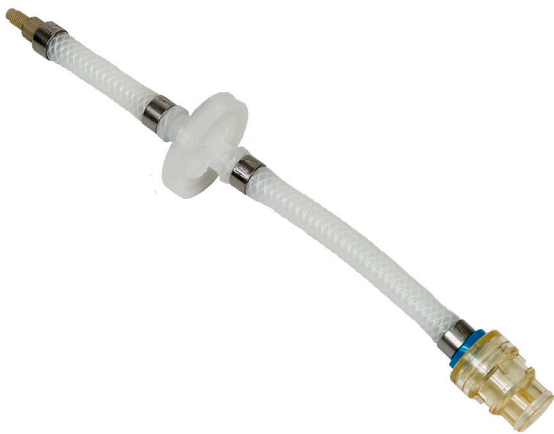

FEATURES

- Disposable
- Quick to remove and replace
- User friendly

SPECIFICATIONS

Measurements:

- Width: 2.4 in. (Max.)
- Height: 10.3 in.
- Depth: 2.4 in. (Max.)
- Weight: 0.11 lbs.

Utilities and/or Components Required:

- Sample Pilot SP100 or SP200